Solis
OF SWITZERLAND

Solis NAIL STUDIO Typ 78

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- 11 = large club shaped milling cutter
- 12 = oval milling cutter
- 13 = pointed grinder
- 14 = conic grinder
- 15 = conical, pointed milling cutter
- 16 = round milling cutter for corn roots
- 17 = long cylindrical wheel cutter
- 18 = buffer with grinding and polishing plate

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Recommendation:

Be sure motor is off when changing the tool.

Remove tool immediately after use.

Clean the milling cutters with a small hard brush from time to time.

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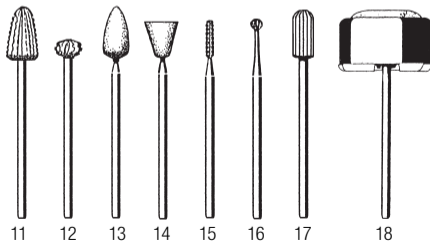
NAIL STUDIO for healthy feet and well cared hands

This instrument is specially developed to give your nails beauty and your skin daily and necessary care carefully and easy thanks to

- its easy manual
- its exclusiv equipment
- its nine NAIL STUDIO precision tools

Recommendation

Study carefully this manual and your NAIL STUDIO will serve you well.

For an improper handling we donot take any responsibility

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Be sure main voltage is identical to the one of the instrument

How to set the tools

- hold the NAIL STUDIO in your hand
- insert the required tool into the opening at the top
- overcome the resistance which occurs at about half the distance and push forward up to the stop


Be sure motor is off when changing the tool!

After use put the tool back to the case

Recommendation

For best results apply on **dry and tight** skin.
Work with **slight pressure** in a counter clock-wise motion.
Use short fast strokes.

Make it a rule to treat your feet periodically.
Do not wait until they hurt.

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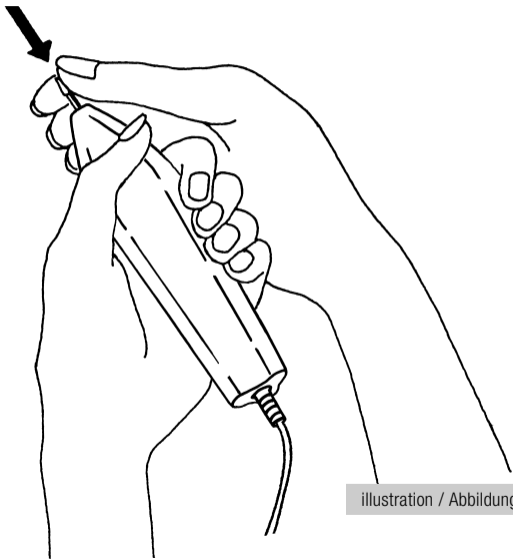How to remove horny skin and callus

illustrations 2 and 3

With your free hand hold the foot to be treated and stretch the skin tight with your thumb.

Tool Nr 12 = **oval milling cutter**
remove horny skin roughly

Tool Nr 11 = **large club shaped milling cutter**
to smooth out afterwards

Such treated the skin will become very smooth

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How to remove corns

illustration 5

Use tool Nr 12 = **oval milling cutter**
Bring foot into a comfortable position.
First remove the edge of the corn root
with soft circling movements
performed without pressure.

illustration 4

The root of the corn (black dot in center)
is removed with special
tool Nr. 16 = **round milling cutter**
by means of slight dabbing movements.

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Treatment of ingrown nails

illustration 6

If the nail has grown into the side bed use tool Nr 15 = **conical, pointed milling cutter**
Cut though the nail sideways and remove the ingrown part with dabbing movements.

Treatment of thickened nails

illustration 7

Use tool Nr 17 = **long cylindrical wheel cutter**
Reduce the surface of the nail to normal thickness with slight fast movements performed without pressure.
You feel it: a normal and thin nail is always elastic and supple.

Recommendation

To avoid inflammations, dab treated parts with oiled cotton.

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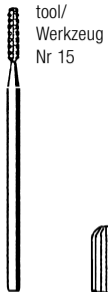
illustration 8

Use tool Nr 14 = **conic grinder**

Hold the instrument as shown on illustration Nr 18

Turn the head of the instrument so that only the nail is touched by the grinding wheel and not the skin. After shaping the nail smoothen the sharp edges.

You notice:

Your fingernail becomes its original form and feels very smooth.

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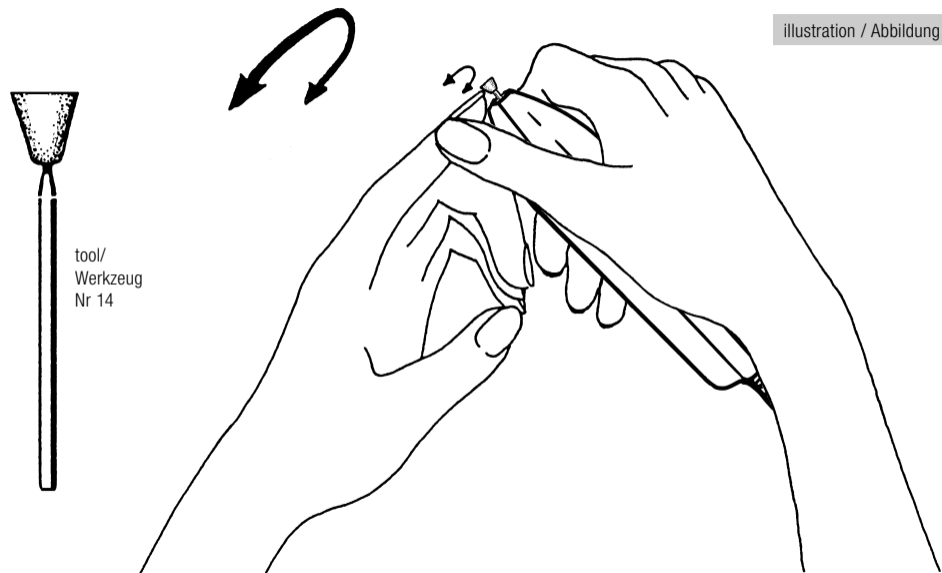
Polishing of fingernails

illustration 9

To polish nails must be clean and dry.
Please remove old nail varnish.
Use tool Nr 18 = **buffer**

step 1

Start to work with the rough grain plate.
Move slightly without pressure over the nail.
Do not remove the nail dust!

step 2

Change to the fine grain polishing plate and
repeat the procedure until the nail shines.
Nails treated this way become strong and healthy
and have a better resistance and a natural shine.

Changing of plates

(cloth holder) and grinding or polishing plate. To change
the plate, pull off the socket. Put the socket upside down
on a flat surface and insert the spare plate.

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chuck /
Hülse

plates /
Schleifplättchen

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The buffer (tool Nr 18)
consists of 3 parts:
chuck, socket and plate

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Removal of cuticle

illustration 10

Remove cuticles only after polishing.

Use tool Nr 13 = **pointed grinder**

Loosen your skin at the root of the nail.

Important: hold instrument so that the surface of the nail is not damaged.

illustration 11

Hard skin on nail sides and under the nails is removed with

tool Nr 15 = **conical, pointed milling cutter**

Recommendation:

Dab treated parts with oiled cotton.

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SOLIS AG